



Super Mini Piston Nebulizer

Model: NI60 Cutie



Valve Adjustable Technology

Our VA technology allows users to adjust different levels of Nebulization rate ranging from 0.08(closed) – 0.2ml/min(fully opened) at consistent particle size less than 2.6 µm.




Mini Pump

The mini pump is a brand new design for small portable nebulizer. The components of the pump are of the highest quality and manufactured by Rossmax Japan. Moreover, the pump is tailor-made which means that the pump's characteristics and functionality can be set and adjusted to meet nebulizer specifications.

Patent Nos :
CN: ZL 2012 2 0443597.3
ZL 2014 3 0053957.3
Other patents pending

- MMAD ≤ 2.6 µm;
Fine Particle Dose (FPD): 75-80%
- Consistent fine particle size for efficient respiratory treatment
- Patented Valve Adjustable nebulizer bottle (VAT)
- Mini pump technology
- Low noise and quiet
- Ideal for kids
- Very compact (215.5 x 134.6 x 93.3 mm) and light weight (334 g)
- Mouthpiece and mask for child included
- 12V adaptor included

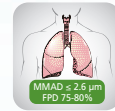

Efficient Treatment



Patented Bottle



Mini Pump



Low Noise Design



Nebulizer Kit, Angled Mouthpiece



Mask for Children



Model	Qty per carton	Carton volume
NI60Q	8 pcs	0.064 cbm/ctn
QTY per 20'	QTY per 40'	QTY per 40'HQ
3496 pcs	7000 pcs	8120 pcs