



Blood Glucose Monitoring System

HS200 Blood Glucose Test Strips



- Blood glucose test strips for self –testing and single use
- For diabetic patients for diabete control
- For quantitatively measuring glucose in whole blood
- 1 vial of 50 strips
- To be used 6 months after opening
- 22 months life span
- Test time: 5 seconds
- Only works with rossmax HS200 glucose meter
- IVD Directive 98/79/EC & EN ISO15197:2015 compliance

Technical Specification

- Measuring method: GDH-FAD Auto coding
- Measuring range: 20~600 mg/dL (1.1 ~ 33.3mmol/L)
- Storage conditions : Temp: 4°C~30°C, do not refrigerate
- Blood source: fresh whole capillary blood
- Volume of blood: min 0.7µl

GDH

GDH-FAD
Technology

FREE


Auto Coding

0.7µl

Measures with
Small Blood Sample

Advanced GDH-FAD Technology

The GDH-FAD test strips can eliminate interference from oxygen variations and non-glucose sugar such as maltose and galactose, providing you with the most reliable and accurate results.

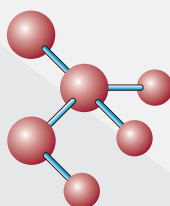

GDH-FAD

Enzyme	GDH-FAD	GDH-PQQ	GOD
Affected by Oxygen	No	No	Yes
Affected by Galactose	No	Yes	No
Affected by Maltose	No	Yes	No