



Blood Glucose Lancing Device



**Painless for
Maximum Comfort**



HIGH DESIGN
QUALITY

reddot design award
winner 2015



DESIGN
AWARD

Taiwan 2017
Golden Pin Winner

- Patented safety switch design
- Patented One-handed lancet ejector
- Compatible with all standard lancets
- Ergonomic design for soft touch
- Easy to use for all users
- 6 Depth setting
- AST (Alternative Sites Testing) cap is available

Ergonomic design enables the comfortable press on finger and lowers the psychological stress before activating it.

Depth Setting: Big window shows the penetration depth with numeric display which is easy for the user to read.


Safety Switch: Patented safety switch design: when the cap is close, it allows only load but lock the function of lancet ejection. When cap is open, it allows lancet ejection but release the trigger automatically to prevent accidental pricking.

Thumb Activator: Big release button, comfortable and easy to operate.

Lancet Ejector: Patented feature of lancet ejection by one thumb pushing with only one hand operation.




Depth adjustable cap

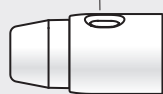

Lancet holder

Release button

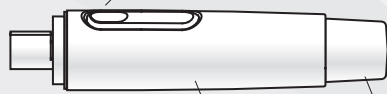

Hub

Plunger



Clear cap (Optional)



The protective cover

The disposable lancet

Model	Qty per carton	Carton volume
LD-E1	540 pcs	0.045 cbm/ctn